

CLINY

“NEW” All Silicone Tracheostomy Tube



Improved Points

1. Cannula/shaft combined two silicone materials

“Harder” rubber which is preventable kinking of the cannula caused by bending cervix.




“Softer” rubber which can reduce the stress to trachea as a measure to prevent formation of granulation.

2. Suction Line

Placed the suction eye at the side of the cannula not to be blocked up the line by contacting to the trachea.



Increased the suction efficiency by the position of the suction eye close to the cuff.

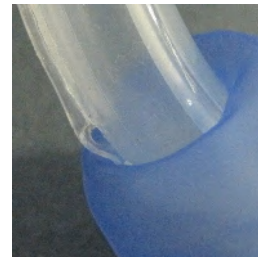

3. Low-leakage Soft Cuff

Lower gas-permeability to maintain the form of the cuff for longer time.
Higher wrinkle-resistance for easier insertion.

4. Multiple sizes

Various sizes are available for selecting appropriate size according to cases.

Specifications

Product Code No.	Size	I.D. (mm)	O.D. (mm)	Length of Cannula (mm)
800 001 0985	22Fr	5.0	7.3	57.0
800 001 0986	26Fr	6.0	8.7	63.0
800 001 0987	30Fr	7.0	10.0	71.0
800 001 1578	32Fr	7.5	10.7	73.0
800 001 0988	33Fr	8.0	11.0	75.0
800 001 1579	35Fr	8.5	11.7	78.0
800 001 0989	37Fr	9.0	12.3	80.0
800 001 0990	40Fr	9.5	13.3	83.0