

# Tierrett Catheter for cats

Registration No. of Medical Equipment for Animal Use: 26-3981

Urinary catheter developed with the concept of "not easily coming out"



- Tube made of pliable polyurethane
- Two types are available: Straight and pig tail
- The pig-tail type outfitted with side holes is suitable for urethral catheterization/placement.









# Tierrett Catheter for cats

## **Features and Advantages**

### Features common to both straight and pig-tail types

- Polyurethane tube
  - The tube is made of pliable and elastic polyurethane.
  - Three X-ray contrast stripes help visualization of the placement under fluoroscopic control.
- Connector
  - A lightweight/compact lure connector is adopted to reduce stress.

## Pig tail

- o Shape of tip
  - The shape of a pig tail provides resistance against the urinary sphincter.
    Thus, it is effective for urethral catheterization on the spot.
  - Loop diameter is set at 10 mm.
- Side holes
  - Three holes (2 round and 1 oval shapes) smoothen urethral catheterization.



# **Product Specifications**

#### ➤Straight (Open tip type)

Product code	Size designation	Outer diameter (mm)	Effective length (mm)	JAN code
800-001-3516	3Fr	1.10	120	4542187191844
800-001-3518	4Fr	1.33	120	4542187191868

5-pack box/Sterilized

#### > Pig tail (Closed tip with side holes)


Product code	Size designation	Outer diameter (mm)	Effective length (mm)	JAN code
800-001-3517	3Fr	1.10	120	4542187191851
800-001-3519	4Fr	1.33	120	4542187191875

3 side holes with stylet; 5-pack box/Sterilized

- ➤ Usage Note
  - When inserting the pig tail, stretch the catheter and insert the stylet.
  - Do not push with force. Stop using the product when having difficulty in inserting it, and use an alternative method.
  - □ Please note that the product specifications are subject to change without notice to reflect the results of our consistent R&D efforts.
  - Read the package insert before using the product.



https://www.createmedic.co.jp/english/

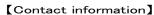

Animal Medical Business Department E-mail: oem-infoch@createmedic.co.jp
[Head Office] 2-5-25 Chigasaki Minami, Tsuzuki-ku, Yokohama, Kanagawa, Japan

